


SARS-CoV-2 CE

Antigen Rapid Test Kit

(Colloidal Gold Immunochromatography)

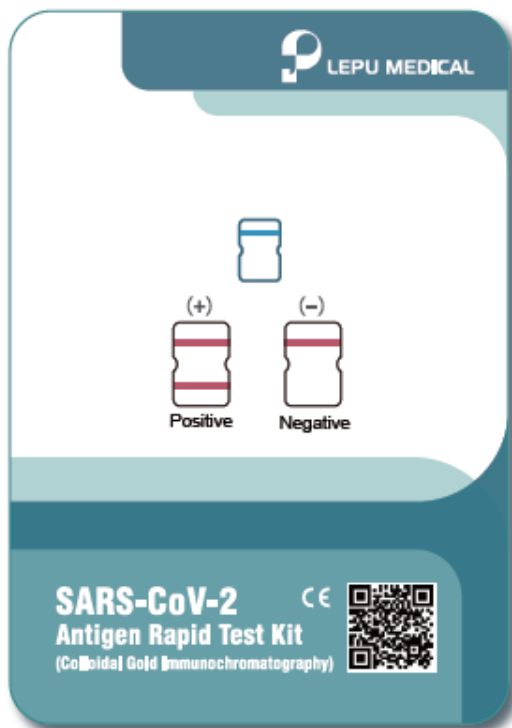

- ① Non-invasive
- ② Simple to use
- ③ Convenient, no devices required
- ④ Rapid, get result in 15 minutes
- ⑤ Stable, with high accuracy
- ⑥ Inexpensive, cost-efficiency



Hospital



Test Site



Airport



Station



Hotel




Corporation



Mass
Screening


Instruction

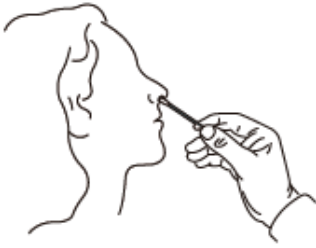


Step 1: Use swab to collect sample.



Step 2: Peel off adhesive.

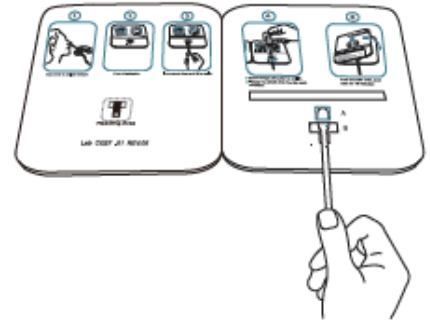

Step 3: Insert swab from well B to well A.

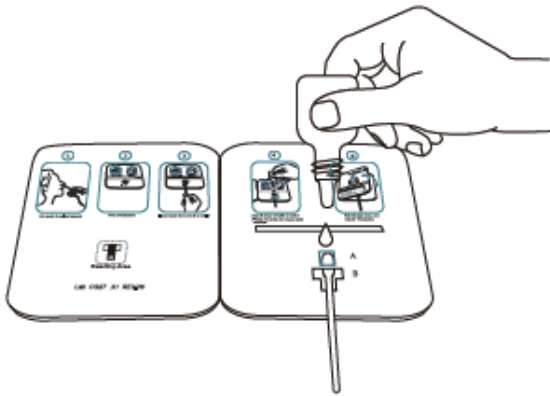

Step 4: a. Add 6 drops of buffer to well A
b. Rotate the shaft, two rounds each direction.



Step 5: Fold left side over, and wait for 15 minutes.

Result Interpretation

Positive

(+)



Negative

(-)



Invalid

